ZINC PHOSPHATE CEMENT

INSTRUCTIONS FOR USE

English
Español
Français

 CONAMCO, S.A. de C.V.
José Ma. Bustillos No. 28-Bis, 1st. and 2nd. Floor.
Col. Algarín, Alcaldía Cuauhtémoc.
P.C. 06880, Mexico City, Mexico.

 Ivonne Richaud.
Travessera de Gracia Street No. 315, 2nd.-1st.
P.C. 08025, Barcelona, Spain.

REV: "0" / 21 JUL 202

ZINC PHOSPHATE CEMENT

ENGLISH

PRODUCT DESCRIPTION: Zinc Phosphate Cement is a material used for cementations of inlays, crowns, bridges and orthodontic appliances and occasionally as a temporary restoration (when it mixes more thicker).

CLINICAL BENEFITS:

- Long history of use
- Excellent compressive strength making the restoration more durable
- Low solubility in the oral médium, which increases its durability.
- Versatil and economic.
- Allows for thin film thickness making it more adaptable and precise.
- Once set, it remains stable.
- Compatible with metal and ceramic.
- It returns the masticatory forces of the treated piece and in cases o orthodontic keeps the devices in position during the treatment.
- Rebuild and return chewing functionality.

PERFORMANCE CHARACTERISTICS:

- Its compressive strength is relatively high and tensile strength is low.
- Adhesion is only mechanical.
- High mechanical resistance.
- Solubilidad moderada.
- Moderate working time.
- Antibacterial effect.

COLOR AND SHADES: Opaque

INTENDED USE For permanent luting cementation of crowns and bridges restorations. (To attach or cement indirect restorations to prepared teeth).

INDICATIONS:

- Permanent cementation of crowns, bridges, onlays and orthodontic appliances.
- Base or lining cement in deep cavities.

- Temporary Restorations if mixed thicker.
- Under amalgam or other metallic restorations for thermal isolation.

INTENDED USER: Dental professionals only with an academic degree, experience, and training for use dental materials, it must be manipulated by a dental professional only.

PATIENT POPULATION: All patients of any age, gender, weight, or ethnicity, with permanent dentition, who, for different reasons, must make restorations with Medental Zinc Phosphate Cement.

PRECAUTIONS:

Some patients can present allergies to the components of the formula, in this case, suspend it use.

Keep away from Sunlight.

Keep away from rain.

Keep away from the reach of children.

Keep away from heat or sparks.

Temperature limit between 4°C and 24°C.

Fragile.

Non-Sterile.

No requires refrigeration.

Only for professional dental use.

Room temperature before use.

Not require cold chains to transport.

Do not use it after the date of expiry.

CONTRAINDICATIONS:


- The material described must not be used when the application techniques are not possible.
- It is not recommended to use as a cavity lining because it can irritate the pulp.
- Contraindicated in: Deep cavities or those close to the pulp because it can cause irritation.
- Patients with pulp hypersensitivity or marked sensitivity.
- It is opaque, so it is not recommended when there is a need for great aesthetics.
- Long-term provisional restorations, as it does not release fluoride.

- Cases where chemical adhesion to the dental structure is required.
- Areas subjected to extreme stress where their fragility may be a disadvantage.

COMPOSITION: Zinc Oxide and phosphoric acid.

INSTRUCTIONS FOR USE:

- Use a cool mixing slab (72°, 22°C) and a non-corrosive spatula for mixing ZINC PHOSPHATE CEMENT.
- Use a ratio of 1.3 to 1.4gr. of POWDER to 0.5cc of LIQUID.
- Divide POWDER into four (4) equal portions. Divide one (1) of these quarter portions into two (2) parts.
- Spatulate 1/8 of the portions into the entire amount of the liquid for 10 seconds.
- Add the other 1/8 of the portions and mix for 10 seconds.
- Add consecutively two (2) of the quarter portions, spatulating 10 seconds after each addition.
- Add the final quarter portion and mix the entire amount for 20 seconds. This will give 1 minute total mixing time.

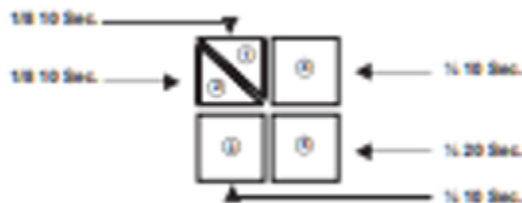MIXING TECHNIQUE:

SETTING-TIME CONTROL:

- Powder-Liquid ratio. Thick mixes set faster than thin mixes.
- The faster the powder is added to the liquid, the quicker the set.

- Prolonged mixing tends to lengthen setting time.

CONSERVATION:

-Close the package immediately after to dispense the product.

-To prevent deterioration of material, keep bottle and jar tightly closed when not in use.

-To eliminate condensation of water vapor on the slab, raise the temperature of the slab slightly if 72°F (22°C) temperature of slab is below the dew point.

-Be certain all the powder is thoroughly incorporated into the liquid in each successive step of the mixing.

WARNINGS and PRECAUTIONS: In case of any change to the original function of the product, stop using it and contact the manufacturer.

-Read carefully and follow all instructions

-May cause mild pulp irritation or inflammation during application.

- It can be difficult to remove if used as temporary cement.

-Before the use of these products always use rubber dam gloves and protective glasses.

-In case of contact with the eye, wash immediately with plenty of water for at least 15 minutes and if the irritation persists consult a doctor.

-In case of contact with skin wash the area with soap and water.

-The product should not be ingested or aspirated. If accidentally swallowing induces vomiting immediately and get medical attention. In case of inhalation, clean your nose with water and seek medical assistance.

-Ask for patient to rinse his mouth immediately after treatment.

WARRANTY:

If the use of this product on a patient or anyone else causes some anomaly, please send detailed report to: Conamco, S.A. de C.V. email: atencionclientes@conamcosacv.com.mx

Warranty: MEDENTAL INTERNATIONAL and/or CONAMCO S.A. de C.V. warranty the product against

any fault of manufacture of the material. MEDENTAL INTERNATIONAL AND/OR CONAMCO S.A. de C.V. does not take responsibility for the use of the product outside of what is specified in this instruction sheet. Being the only obligation is the replacement of defective product.

The summary of Safety and clinical Performance (SSCP) can be retrieved from the European Database on Medical Devices (EUDAMED) at:

<https://ec.europa.eu/tools/eudamed>

DISPOSAL CONSIDERATIONS:

New, clean nitrile gloves should be worn.

1) *Removal of excess.*

- Fresh (before complete setting): Remove with a probe or spatula. recommends cleaning immediately.

- Setting: Using cement chisels, curettes or metal excavators with controlled movements, pass waxed dental floss in interproximal areas, pulling downwards to avoid fracturing restorations.

2) *Removal under restorations or crowns.*

- It is not possible to remove it without removing the restoration; in these cases, the crown must be sectioned with a carbide/diamond bur and the tooth cleaned manually or with ultrasound.

*In both cases, it is recommended to protect soft tissues, work with good lighting, and use a magnifying glass to identify cement residues.

3) Dispose of with solid waste in accordance with the regulations of each country.

FEATURES/CLAIMS: Water Base Cements are a material used for cementations of inlays, crowns, bridges and orthodontic appliances and occasionally as a temporary restoration. It is prepared by mixing Zinc Oxide and magnesium oxide powders with a liquid consisting principally of phosphoric acid, water and buffers. This Cements can be used for fixing crowns, bridges, posts, cores, in various materials with clinical properties like modulus resilience, mechanical adhesive, solubility, adhesion and retention, long term survival event rates. Water Base cements are effective

in clinical context and in short-term clinical use, and they can increase strength with storage.

Zinc phosphate cement is recommended where maximum mechanical retention is required, is superior for retention, tensile bond strength, but it can reduce strength with humidity. To take good care of the restoration, it is recommended to avoid biting on very hard things so that this material does not splinter. Also moderate the consumption of highly pigmented foods such as coffee and wine, as well as maintain good dental hygiene and use dental protector if the patient has Bruxism.

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